

NEW SYRINGE ASSEMBLY, PLEASE READ INSTRUCTIONS BEFORE PROCEEDING

CONTENTS:

Single-use glass syringe prefilled with 0.3 mg of Macugen®.*
Sterile packaged BD® single use 30 gauge x ½" Precision Glide® Luer Lok® needle.

Note: Prior to injection and to ensure proper dose, plunger stopper must be properly aligned so that the top edge of the 3rd rib on the plunger stopper is aligned with the pre-printed black dosing line. See Fig. 2 below for details.

ASSEMBLY:

The injection procedure should be carried out under controlled aseptic conditions, which includes the use of sterile gloves, a sterile drape, and a sterile eyelid speculum (or equivalent). When ready to assemble syringe and administer injection, carefully peel open pouches, remove contents, and place on sterile field. If upon opening the pouch, if the plastic clip is missing or not attached to the syringe, do not use the syringe.

To avoid compromising the sterility of the product, do not pull back on the plunger rod.



- 1 Remove the syringe from the plastic clip.
- 2 Twist off cap.
- 3 Attach the sterile administration needle (included) to the syringe by screwing it into the syringe tip. —Another sterile administration needle may be used in lieu of the one included. Remove the plastic needle shield from the needle.
- 4 Holding the syringe with the needle pointing up, check the syringe for bubbles. If there are bubbles, gently tap the syringe with your finger until the bubbles rise to the top of the syringe. SLOWLY depress the plunger rod to eliminate all the bubbles and to expel the excess drug so that the top edge of the 3rd rib on the plunger aligns with the pre-printed black dosing line (see Figure 2, right).
- 5 Inject the entire contents of the syringe.

PRIOR TO INJECTION


Fig 1. Before expelling air bubble and excess drug



(Actual air bubble formation may vary)

READY FOR INJECTION

Fig 2. After expelling air bubble and excess drug

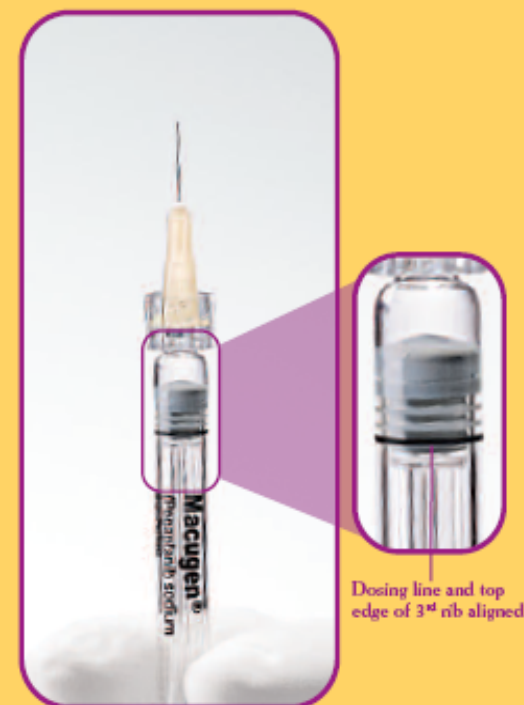

Dosing line and top edge of 3rd rib aligned

*Each 1-mL syringe delivers a unit dose of 0.3 mg pegaptanib sodium in 90 µL when administered as directed.

Please see accompanying package insert for full prescribing information.

Syringe is single-use only. Do not reuse or refill. After use, dispose all syringe parts according to standard biohazard disposal procedures.

BD and Precision Glide Luer Lok® are registered trademarks of Becton Dickinson & Co, Franklin Lakes, New Jersey 07417